

# Wobble-not™


## Low Insertion Force Serological Pipet

See the difference. Feel the difference.

### TRY IT

### COMPARE IT

- Compare the force you need to insert
- Compare the amount of wobble
- Compare the amount of dripping

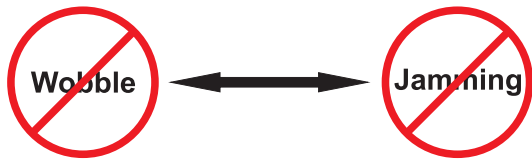

**Low Insertion Force**

Less Wrist Strain  
 Less Pain

**More Ergonomics!!!**

All pipette controllers have the same basic nozzle design for holding onto serological pipets.



#### Wobble-not™



With low insertion force, the Wobble-not pipet is held stably



#### Standard Pipet

If you don't use enough force, the plug end is not held tightly enough and the pipet is more wobbly.

To get the pipet to wobble less, you must use much more force to insert the pipet.

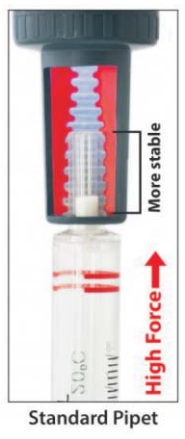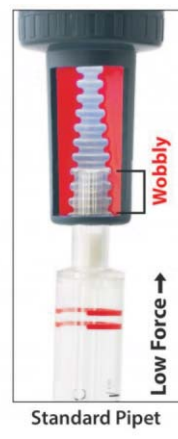


**Wobble-not = less pain and maximum stability**

# Low insertion force = Less wrist strain = Less pain = More ERGONOMIC!

Wobble-not serological pipets are a revolutionary, new design and the only serological pipets that will give you two zones of contact with relatively low insertion force. No jamming! And furthermore, this results in maximum stability for your pipetting and LESS DRIP! No other serological pipet comes close to achieving the same amount of Wobble-not!

## Less Sag = Less Wobble = Less Drip!

### Why it works

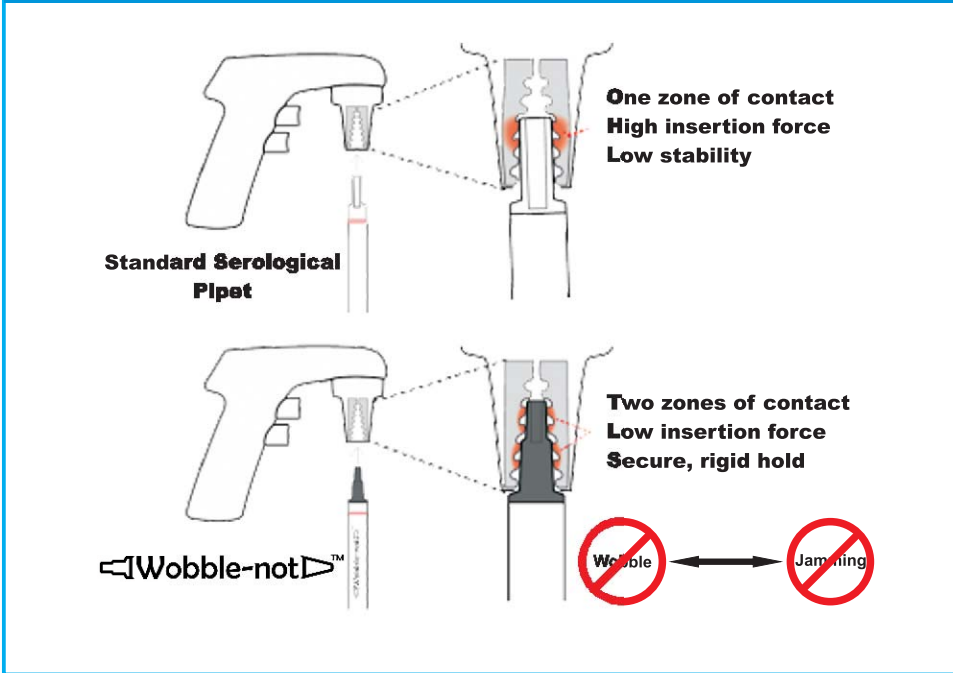

Wobble; wob·ble (w? b? l); noun. A movement or rotation with an uneven or rocking motion or an unsteady motion from side to side.

Traditional serological pipets have a straight cylindrical shape vs. the Wobble-not, which has a two-tiered step design. This helps to decrease insertion force, while increasing pipet stability in your controller. So, no jamming and less dripping!

### The Difference Is Crystal Clear!

There Is No Comparison. Wobble-not is superior in ergonomics and performance compared to other standard serological pipets.

- » Unique patented 2-tier step design at plug end
- » Individually wrapped in paper/plastic
- » Sterile – RNase/DNase free, Non-pyrogenic
- » Accuracy  $\pm 2\%$  at full volume
- » Filter at plug end made of porous polypropylene
- » Ultra-clear markings

Product Code	Volume	Gaduation Interval	Negative Graduation	Quantity Per Pack	Quantity Per Case
4090-0005	5ml	1/10ml	3ml	50pcs.	200pcs.
4090-0010	10ml	1/10ml	3ml	50pcs.	200pcs.
4090-0025	25ml	2/10ml	8ml	25pcs.	200pcs.
4090-0050	50ml	1/2ml	8ml	25pcs.	100pcs.

**Catalyst Biotech**

C-302, Third Floor, Anupam CHS Ltd.  
Plot No.: 79, Panvel - 410 206.  
Maharashtra - India.

E-mail : [catalyst@catalystbiotech.com](mailto:catalyst@catalystbiotech.com)  
Website : [www.catalystbiotech.com](http://www.catalystbiotech.com)

Our Authorised Distributor :

